


This gel depicts the NEB 1kb ladder and a digestion of the plasmid pMYO127TAG-SepT with the following enzymes: **NotI-HF**, **BamHI-HF** and **NotI-HF/BamHI-HF** double digest

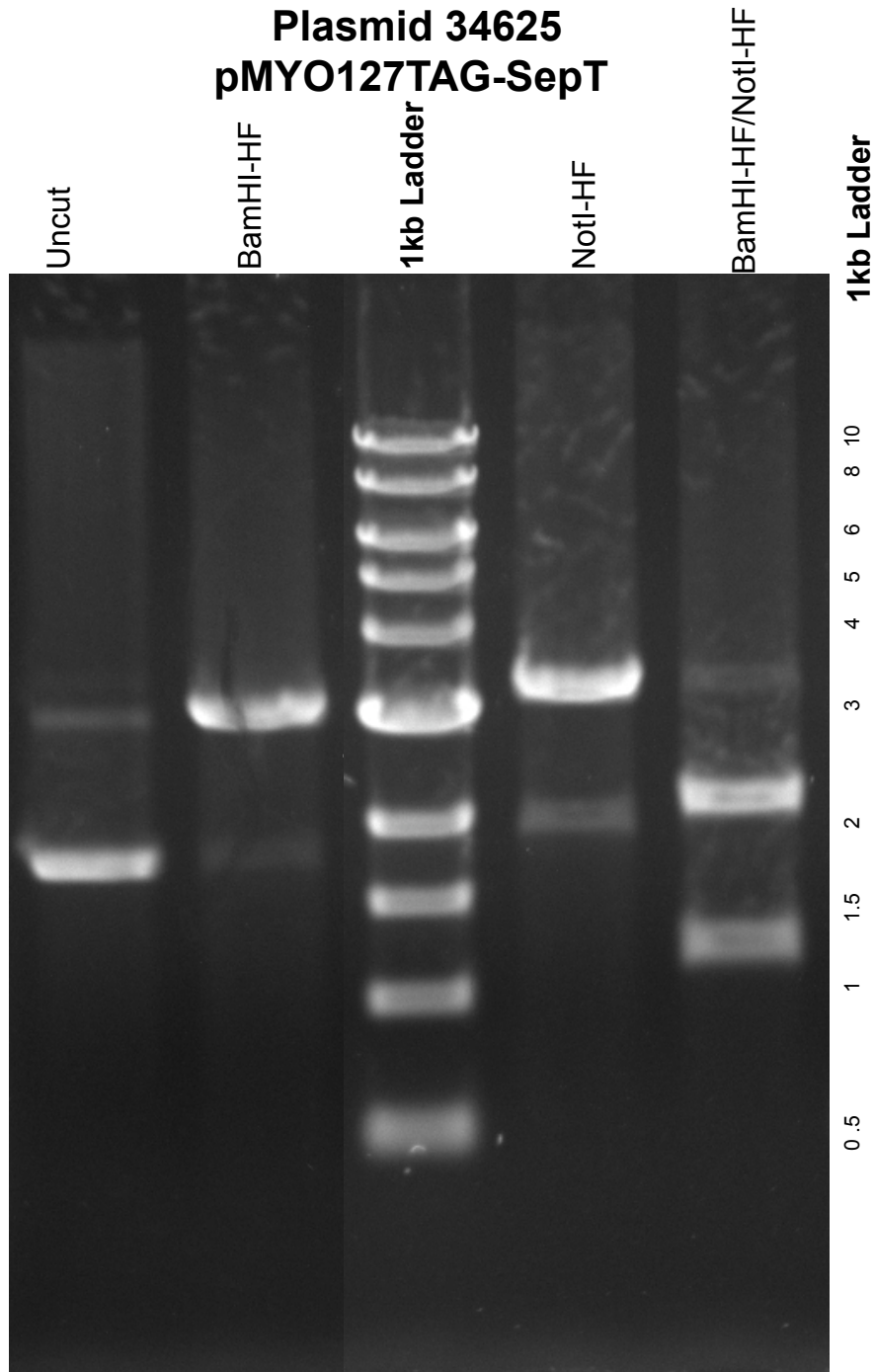

DNA digested with **NotI-HF**, **BamHI-HF** and **NotI-HF/BamHI-HF** double digest for 1 hour

Gel: 0.8% TAE stained with 100 ml 1X TAE and 5uL of ethidium bromide

Ran for 1 hour and 20 minutes at 80 volts